


Connecting Options for Analytical HPLC Columns, Short Columns and Guard Columns with 2 - 3 - 4 - 4.6 mm ID

1. Guard Column and Analytical Column

- 1.1 Connection via capillary or one-piece PEEK column connector (P/N 7.590253)
 1.2 Direct connection via coupling nut (P/N 6.424204) and stainless steel connector (P/N 6.424404), details refer to catalogue 2020-d, page 122.



7.590253

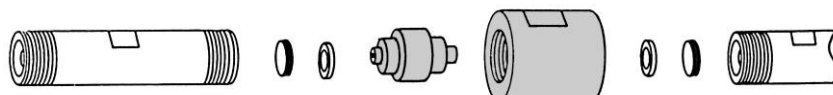


6.424204 + 6.424404

2. Column Connector for ECO Columns

The ECO column is delivered without column heads and fitting adapters, only with end caps.
 For column connection you need:

2 pcs. column head, P/N 6.424201


2 pcs. fitting adapter, P/N 6.424260

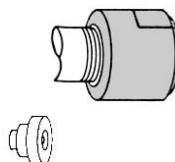

For details refer to catalogue 2020-d, page 122

3. Guard Column Holder and Analytical Column (details refer to catalogue page 123)


- 3.1 Guard column holder with UNF 10/32 female (connection via capillary tubing), P/N 6.424271 and 6.424272

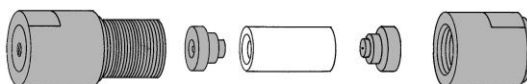

- 3.2 Guard column head, direct connection of guard column with analytical column, P/N 6.424262, 6.424266, and 6.424267



Further info and details can be found in our catalogue 2020-d, page 122-123.

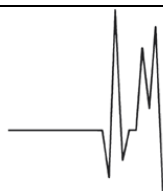

ZIEMER Chromatographie

Am Parir 27 · 52379 Langerwehe / GERMANY
 Phone +49 24 23 / 40493-5 · Fax +49 24 23 / 40493-88
 e-mail: info@ziemer-chromatography.com

Visit our **Online-Shop** under: www.ziemer-chromatography.com